

# Breathing for the most critical patients


**Contact**

**AIR LIQUIDE MEDICAL SYSTEMS**  
Parc de Haute Technologie  
6, rue Georges Besse  
92182 Antony Cedex

Phone: +33 (0)140 96 66 00

Customer service:

**+33 179 51 7001**

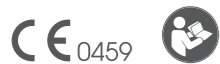

[www.device.airliquidehealthcare.com](http://www.device.airliquidehealthcare.com)

Monnal™ T60 Advanced is a class IIb air turbine medical device for the treatment of infants (from 3 kg), children and adults, intended for hospital use in intensive care, emergency and intra/extra-hospital transport - The acts carried out with Monnal™ T60 Advanced are covered by the health insurance bodies in certain situations.

Manufactured by Air Liquide Medical Systems - Conformity assessment: GMED 0459 - Commercial code: R.C.S. Nanterre 348 921 735 - Please read the user manual carefully.

Air Liquide Medical Systems is committed to an environmental approach with a quality-environment management system certified ISO 14001:2015. For more information, go and visit our website <https://www.device.airliquidehealthcare.com/our-commitments>



VNT04412 - May 2021 version - Non-contractual document - product may change at any time - Photos: ©ALMS, ©Anthony Guerin - Design: CHAM1TH



[www.device.airliquidehealthcare.com](http://www.device.airliquidehealthcare.com)

## Designed for the most critical patients



Monnal™ T60 Advanced is the solution designed to meet caregivers needs and most critical patient management.

Monnal™ T60 Advanced is designed to provide invasive or non-invasive ventilatory support to all patients with acute respiratory distress, and is adaptable to specific environments and situations.

Monnal™ T60 Advanced adapts to all **intra-hospital** environments (from bedside to patient transfers):

- in the intensive care, post-resuscitation and post-operative monitoring departments
- in emergency departments, in the emergency resuscitation room
- during Inter-hospital transfers

Monnal™ T60 Advanced is also suitable for **out-of-hospital transport**, meeting the requirements of emergency situations:

- in land and air mobile units, helicopters, airplanes, and long distance medical transport



## Functions adapted to each department


### Integration of the APRV mode

Monnal™ T60 Advanced mode integrates the **APRV** (Airway Pressure Release Ventilation) mode, enabling invasive ventilation at two asynchronous pressure level

Allows the patient to breathe spontaneously in a pressure-controlled mode

Facilitates spontaneous ventilation of patients while reducing the use of curares<sup>1</sup>

Prevents alveolar distension by spontaneous ventilation and adjustment of the two pressure levels<sup>1</sup>

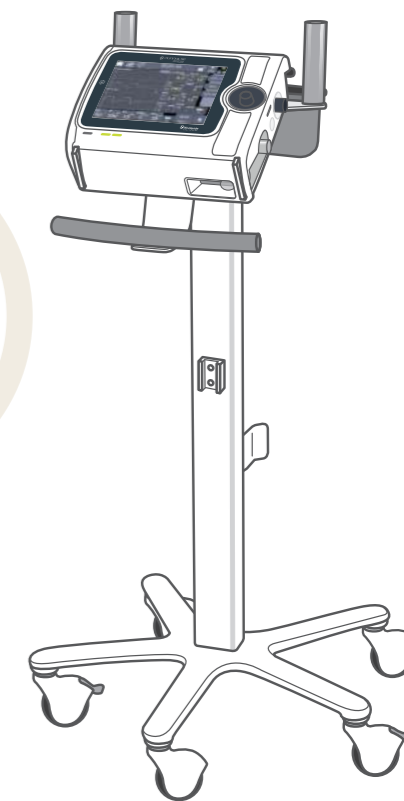

1. Habashi. Other approaches to open-lung ventilation: airway pressure release ventilation 10.1097/01.ccm.0000155920.11893.37

### Advanced monitoring: the partner of secure ventilation

Monnal™ T60 Advanced provides **comprehensive and advanced monitoring**, which is crucial for monitoring and adapting the care of the most critical patients. Thanks to its monitoring functions, the **Monnal™ T60 Advanced** provides:

#### Protective ventilation monitoring

- **VT<sub>e</sub>/PBW:** monitoring of tidal volume exhaled according to the patient's predicted weight (gender and height entered)
- **P<sub>plateau</sub>:** measurement of plateau pressure with alarm thresholds
- **PEEP:** measurement of total PEEP, allowing the estimation of autoPEEP
- **Driving Pressure:** monitoring of the patient's elastic pressure (P<sub>plateau</sub> – PEEP)
- **Static and dynamic compliances:** measurement of static compliance (C<sub>stat</sub>) and estimation of dynamic compliance (DynC)

#### Monitoring of spontaneous ventilation


- **Expiratory minute volume (Spont MVE):** evaluation of the patient's spontaneous breathing level.
- **Spontaneous respiratory rate (Spont RR):** assessment of the patient's spontaneous breathing rate level

# The companion of health care teams

## ⊙ A simple and intuitive interface

**Monnal™ T60 Advanced** is an **intuitive ventilator** that allows teams to quickly become familiar with the ventilator and focus on other aspects of patient management.

The **customizable interface** allows to adapt according to each department's protocols.

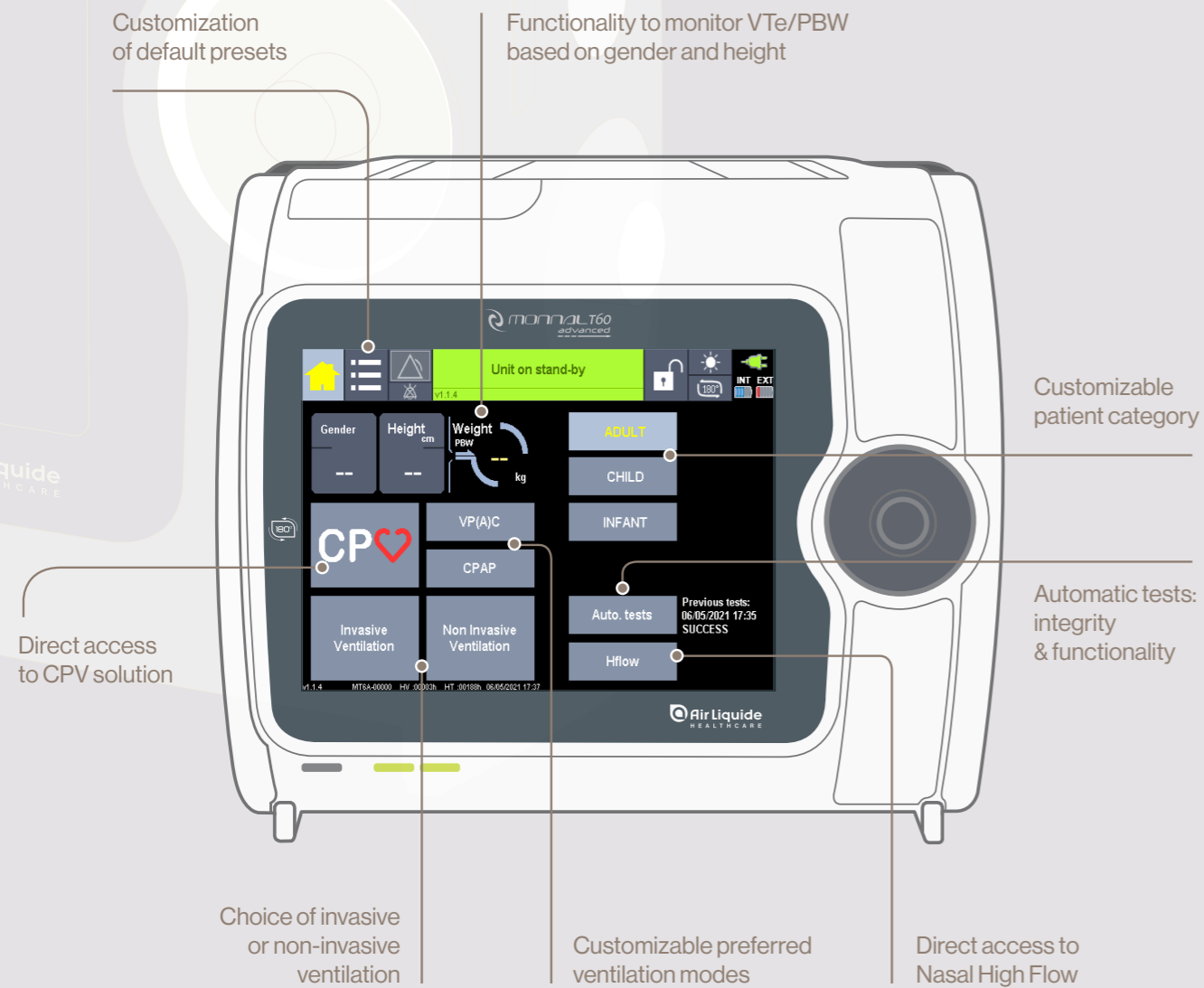

## ⊙ Integration of the CPV solution to manage cardiopulmonary arrest

The CPV solution offers technical innovations for each stage of cardiopulmonary resuscitation (CPR) that are specifically designed for rescuers<sup>2</sup>. It is intended for mobile emergency units, intra- or inter-hospital transport, and specialized departments.



The CPV solution with EtCO<sub>2</sub> (optional) allows monitoring of the patient's capnography while deducing a prognostic factor of spontaneous return of circulation<sup>3</sup>.



## ⊙ ECLS environment

The ECLS (ExtraCorporeal Life Support) function is a **dedicated environment** that supports the clinician in changing ventilation procedures while setting up extracorporeal circulation.



2. D Perkins. European Resuscitation Council Guidelines for Resuscitation 2015 10.1016/j.resuscitation.2015.07.015  
 3. Touma et al. The prognostic value of end tidal carbon dioxide during cardiac arrest 10.1016/j.resuscitation.2013.07.011



## Mobile, autonomous and versatile

### ⊙ Easily transportable

- **Light and compact:** 4 kg
- **Self-powered:** 5 hours in total (internal battery + removable battery, in standard configuration on adult ventilation). The internal battery is only a backup power supply.
- **Built-in accessories:** fixing and protection systems for optimal patient configuration and movement: wall-mounted charging station, carrying bag, universal intra-hospital support, rolling stand

### ⊙ Self-contained air device

- **Turbine technology**
- **Peak flow:** 230 l/min for effective NIV breathing assistance
- **Rapid pressurization (max of 120 cmH<sub>2</sub>O/s):** to best meet the patient's inspiratory demand
- **Quiet operation:** to preserve the noise environment (48 dBA at 1 m)

## Health crisis and extreme situations

### ⊙ Beyond department walls

In a congested ICU situation, **Monnal™ T60 Advanced** extends allows ventilatory management beyond the ICU and recovery room alone.


**Monnal™ T60 Advanced** does not require a compressed air supply. Only (high or low pressure) oxygen, to enrich the mixture.

Thanks to its electrical autonomy and turbine, ventilation is maintained in case of a loss of oxygen.

### ⊙ Nasal high flow oxygen therapy

**Monnal™ T60 Advanced** incorporates the Hflow function that allows caregivers to perform high flow nasal oxygen therapy.

The Hflow function is available for **all patient** categories directly on the **Monnal™ T60 Advanced** start screen.

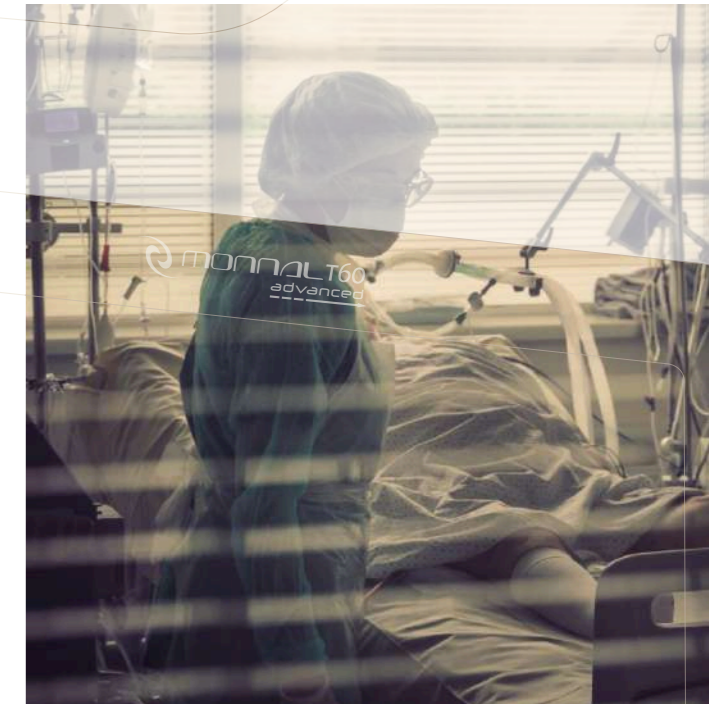

## Monnal™ T60 Advanced For the most critical patients, in extreme situations

**For most critical ill patients**  
ARDS, COVID+, ECMO...

**Specific environments**  
Medical evacuation  
Inter-hospital transport  
In land and air mobile units  
(train, airplane, helicopter)

**Extreme situations**  
Emergencies  
Health crisis

