

JULY 2023









LIGHT DRIVE 2.1


Designed for all manual wheelchair users, whether occasional, temporary or permanent, LIGHT DRIVE 2.1 turns manual wheelchair into ultralight and high performance powerchair in seconds.


Benoit Systemes has been fitting power add-ons on manual wheelchairs in France since 1994. Light Drive 2.1 offers you the possibility to effortlessly add power to your manual wheelchair.



HOW

DOES IT WORK?

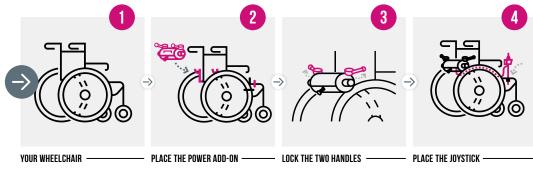

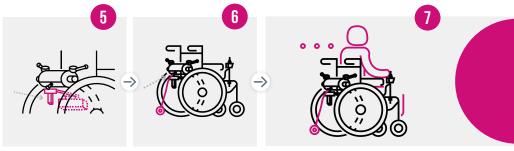

SLIDE THE BATTERY — — INSERT THE ANTI-TIP LET'S GO!




Disengage the power add-on in one movement and feel free to selfpropel.

Jean-Claude

User since 2016 - diagnosed

with multiple sclerosis in 2013.

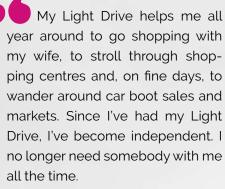





ADVANTAGES OF

LIGHT DRIVE 2.1









INDOORS/

OUTDOORS



LIGHTWEIGHT AND EASILY TRANSPORTABLE

ADAPTS TO ALL WHEELCHAIRS*

INSTANT INSTALLATION

RANGE UP TO 20 KM

*see the compatibility guide



CAN BE **TRANSPORTED** BY PLANE*

Optional 11.4 Ah lithium-ion battery



MADE **REMOVABLE** IN FRANCE



COLOUR OF YOUR CHOICE

*Optional 11.4 Ah battery complies with IATA regulation.



KEY FEATURES



Ergonomic clutch system



Retractable joystick



Anti-tip for safety



Extended colour range



Battery with charge indicator

WIDE RANGE OF COLOUR



CLASSICS

- BLACK
- WHITE
- BLUE
- RED
- KHAKI

METALLICS

- MIDNIGHT BLUE
- GREY
- SILVER
- COPPER
- CHAMPAGNE **GOLD**

TRENDY

- SKY BLUE
- PINK
- GREEN
- YELLOW
- ORANGE
- BLACKCURRANT **OF BURGUNDY**

LIGHT DRIVE MINI 2.1



A smaller version for children

Light Drive Mini 2.1 can be refitted as the child grows allowing to change wheelchair without changing the power add-on.

Wheelchair equipped with a Light Drive maintains its original manoeuvrability and foldability, it is also easy to transport in any vehicle.

Thanks to Light Drive Mini 2.1, your child will be able to ride easily on all kinds of terrain (steep slopes, grass, gravel).

Reverse gear makes the pedestrian pavement and small step more accessible.



TECHNICAL

DETAILS


Our power add-ons can be fitted on most of wheelchairs starting from 32cm seat width and 40 cm between the wheels.

	Maximum user weight	Power add-on weight
Light Drive 2.1	120 kg	8,8 kg
Light Drive Mini 2.1	120 kg	8,8 kg
Light Drive Plus 2.1	240 kg	16 kg

Max speed ———————	—— 10 km/h programmable
Motor —	2 x 130W motors with brakes
Motor Light Drive Plus 2.1 —————	4 x 130W motors with brakes
Climbing ability ————————————————————————————————————	6° (11%)
Manufacturer's warranty	24 months for the power add-on
•	and the battery



LIGHT DRIVE PLUS 2.1

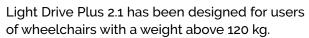

Equipped with four 130 W motors, it will allow you to enjoy the advantages of lightweight power add-ons.

Light Drive 2.1 Plus can be fitted on all foldable wheelchairs. It is provided with a battery of 18Ah and offers a range up to 20 km.



Compatible wheelchairs: a compatibility guide is available to know the models of wheelchairs on which Light Drive 2.1 can be installed, in compliance with general requirement for safety and performance stated in the European Regulation related to medical devices. Find out more





